

Binkie HT **Manual**




1. Introduction

- Product Name: BinkieHT
- Manufacturer: PAPRICALAB Co., Ltd.
- Purpose:
BinkieHT is a heating device specifically designed for BinkieRT tongue positioners to optimize heat transfer for shaping.

2. Parts Description



[Component]

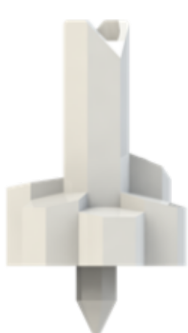

[Holder]



[Usage Example]

- Name: Holder
- Purpose:
Adjust the position of BinkieRT® to optimize heat transfer to the Jaw Positioner.
- Made of highly impact- and heat-resistant materials for long-term use.

3. Function

	Display Panel - [READY] - [HEATING] - [HEATING DONE]
+10 -10	Time Setting Buttons - Increase 10 seconds - Decrease 10 seconds
	Jaw Positioner Heating Button - Preset 30 seconds
	Tongue Positioner Heating Button - Preset 60 seconds
	Pause/Stop Button - Single press: Pause - Double Press: Stop * After pausing, press the heating button to resume

4. Operating Procedure

1. Turn on the power and verify that the [READY] message appears on the display panel.
2. For heating the Jaw Positioner: Attach the holder to the Tongue Positioner before starting the heating process.
3. For heating the Tongue Positioner: Attach it directly to the rotating disc before proceeding.
4. Confirm that the device is operating while the [HEATING] message is displayed.
5. After heating is complete, verify the [HEATING DONE] message and the alarm sound, then remove the BinkieRT®.

5. Specifications & Features

- This product is specifically designed as a heater for BinkieRT®.
- The outer casing is constructed from sheet metal (steel) to provide enhanced durability and protection.
- Dimensions: 350 mm (L) × 120 mm (W) × 250 mm (D)
- Weight: 5.4 kg
- Rated power: 220V / 2000W
- The main unit and internal components are not water-resistant.

6. Safety Precautions

- This product is not water-resistant. Prevent moisture from entering the unit.
- Do not shake the product or subject it to strong impacts.
- Do not disassemble, repair, or modify the product without authorization. Doing so may result in electric shock or malfunction.
- The hot air outlet becomes hot during and immediately after operation. Avoid contact to prevent burns.
- Operate the product only with the specified voltage rating.

 **BinkieHT**

 **Paprica** LAB
Value up Together

PAPRICALAB Co., Ltd.

Headquarters | 2nd floor 11, Yulgok-ro 14-gil, Jongno-gu, Seoul, Republic of Korea

Research Institute | 3rd floor 11, Yulgok-ro 14-gil, Jongno-gu, Seoul, Republic of Korea

Website | www.papricallab.com

E-mail | papricallab@papricallab.com

Phone | +82 (0)2-742-6161